


# Temperature Protocol



GR28008018 28°C to 8°C (18 hour)

Temperature Protocol

	Start Temp °C	End Temp °C	Start Time Min	End Time Min	Description	Time Hr
Step1	-	28	0	30	Prep-Time	0.5
Step2	28	20	30	120	Ramp-down	2.0
Step3	20	15	120	480	Ramp-down	8.0
Step4	15	8	480	1080	Ramp-down	18.0
Step5						

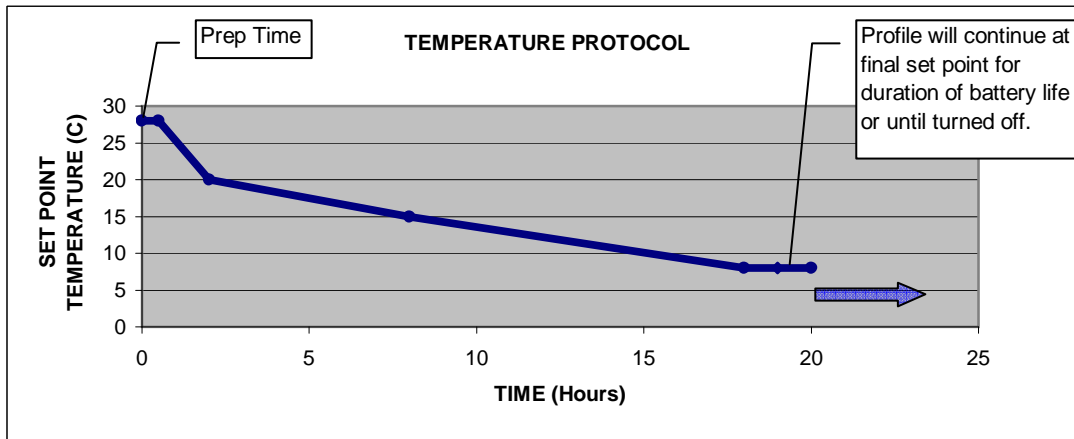

www.microQ.com

Temperature Protocol for Transport of Chilled Oocyte & Embryo

Temperature Protocol for Transport of Chilled Semen

Genesis Series Shipping Devices & Smart Chips