


# Temperature Protocol



## GR28015018 28°C to 15°C (18 hour)

Temperature Protocol

	Start Temp °C	End Temp °C	Start Time Min	End Time Min	Description	Time Hr
Step1	-	28	0	30	Prep-Time	0.5
Step2	28	22	30	120	Ramp-down	2.0
Step3	22	18	120	480	Ramp-down	8.0
Step4	18	15	480	1080	Ramp-down	18.0
Step5						

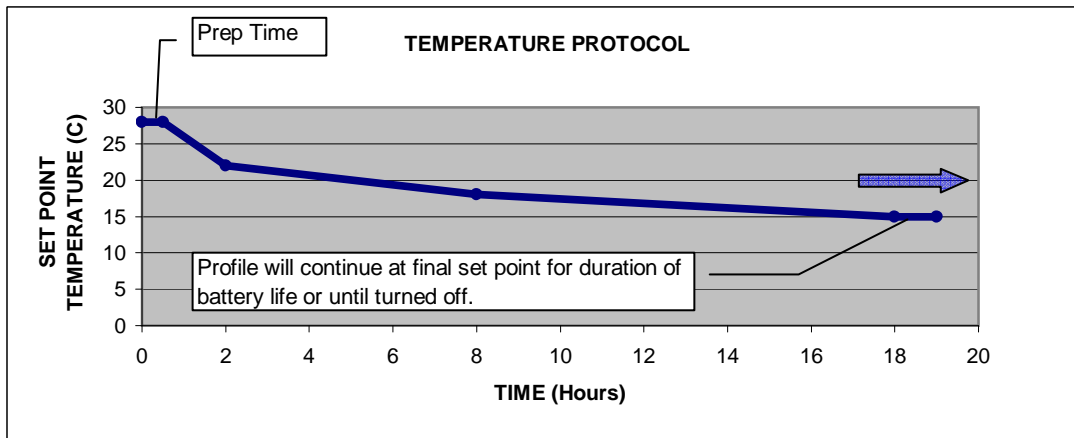

www.microQ.com

Temperature Protocol for Transport of Chilled Oocyte & Embryo

Genesis Series Shipping Devices & Smart Chips